

Passy Muir Valves

All Passy Muir[®] Valves

- Have the patented closed position “no leak” design
- Are designed to fit the universal 15mm hub of disposable tracheostomy tubes including neonatal and pediatric tubes
- Can be used interchangeably by tracheostomized and ventilator dependent patients with the exception of the PMV[®] 2020



PMV[®] 007
(Aqua Color™)
Designed for ventilator application with disposable tubing



PMV[®] 2000
(clear)
Low profile design with lower resistance, recommended for outpatient use



PMV[®] 2001
(Purple Color™)
Low profile design with lower resistance, recommended for inpatient use



PMV[®] 005
(white)
David's Original Design

Passy Muir Valve Accessories

Passy Muir[®] Valve Oxygen Adapter PMA[®] 2000

The PMA[®] 2000 Oxygen Adapter snaps onto the PMV[®] 2000 and PMV[®] 2001 and allows for easy delivery of low flow supplemental oxygen and humidity.



PMV[®] 2001 with PMA[®] 2000


PMV[®] 2001
(Purple Color™)

PMA[®] 2000



PMV[®] Secure-It[®]

The Secure-It[®] strap is designed for use with the lower resistance 2000 series valves and attaches the valve to the trach tie to prevent loss.

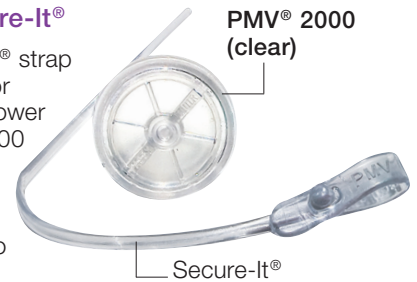


Passy Muir Valves for Metal Tubes

PMV[®] 2020 (clear)

The PMV[®] 2020 (clear) with the PMA[®] 2020-S adapter are for use with metal Jackson Improved Tracheostomy Tubes (Sizes 4, 5, and 6).

PMV[®] 2020
(clear)

PMA[®] 2020-S

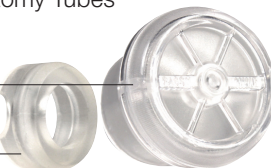

PMA[®] 2020-S




PMV[®] 2020
(clear)

Jackson Improved Tracheostomy Tube

PMV[®] 2001
(Purple Color™)

The Jackson Improved Tracheostomy Tube with Permanent 15mm Adapter can accommodate all Passy Muir[®] Valves except the PMV[®] 2020.

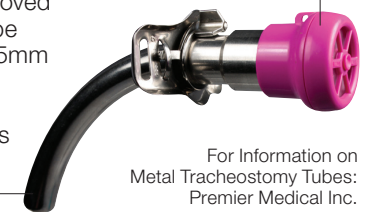

For Information on Metal Tracheostomy Tubes:
Premier Medical Inc.
www.premusa.com
888-670-6100

Passy Muir Adapters

Ventilator Adapters

Passy Muir Adapter

Part # PMV-AD1522
15mm x 22mm Adapter



Passy Muir Adapter

Part # PMV-AD22
22mm Silicone Adapter



Adapters for Duckbill Style Tracheostomy Tubes*

DB15™ Adapter

Part # PM-DB15
15mm Adapter









DigiSil™ Adapter

Part # PM-OCLN
15 mm Silicone Adapter
Includes Secure-It[®]






* Not for use with mechanical ventilation

Ventilator Connections

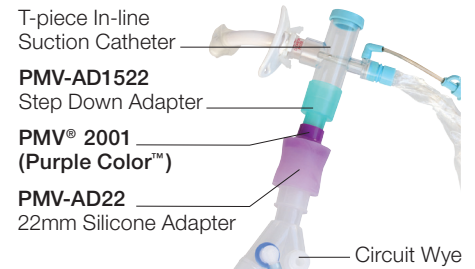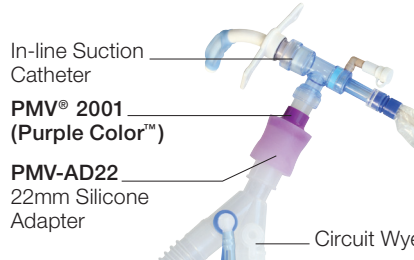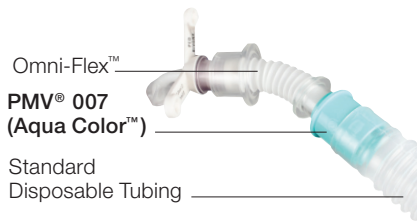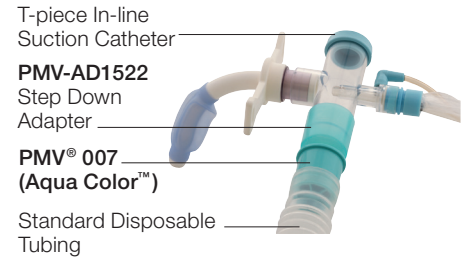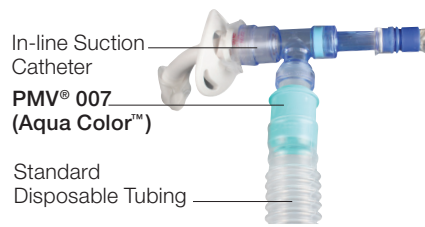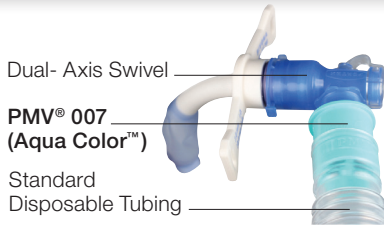


Pediatric Ventilator Connections

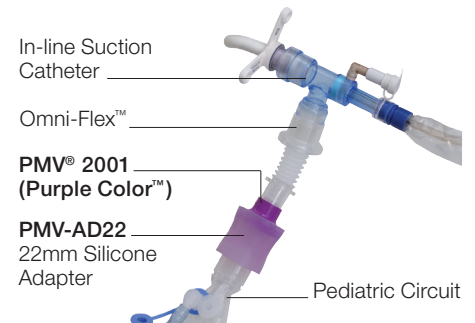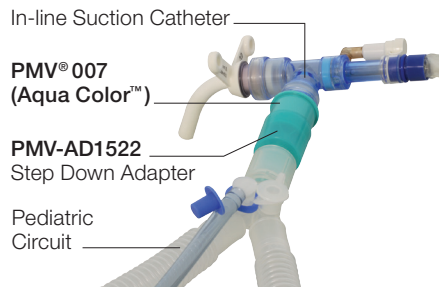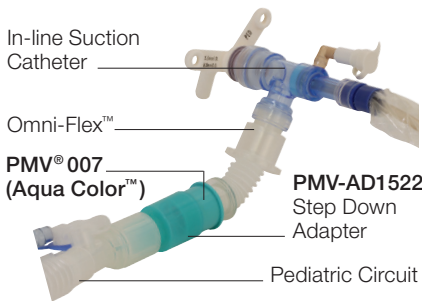

Care & Cleaning of the PMV[®]

Daily Cleaning Procedures

1. Swish valve in pure soap and warm water.
2. Rinse valve thoroughly in warm running water.
3. Allow valve to air dry thoroughly before placing it in storage container.
4. DO NOT use hot water, peroxide, bleach, vinegar, alcohol, brushes, or cotton swabs to clean valve.

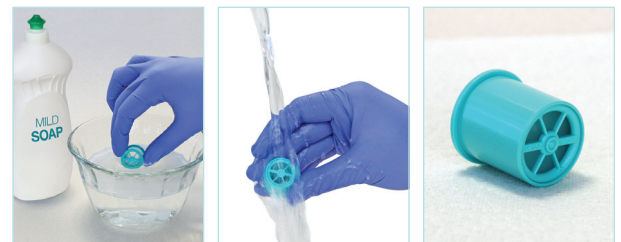

Education & Clinical Support



Speech and Respiratory Clinical Specialists are available to answer your questions

1.800.634.5397



Email your questions

info@passymuir.com

Educational Presentations



- Facility Locations
- State and National Conferences
- Society Meetings
- Seminars
- Colleges and Universities



- Live Special Events Webinars
- Self-Study Webinars
- Free CEU Courses
- Free Clinical Resources

www.passymuir.com

PASSY MUIR PRODUCTS
PROUDLY MADE IN USA

